

ESCOBAR - URIBE ASOCIADOS

A B O G A D O S

Carrera 11A No. 94 - 76 - Of. 305 - P.O. Box 94948, Bogotá, D.C. - Colombia

Tels: (571) 6162051 - Fax (571) 6161889

www.escobaruribe.com - info@escobaruribe.com

53/

2956/P.

SUPERINTENDENCIA DE INDUSTRIA Y COMERCIO DIVISION DE NUEVAS CREACIONES

E. S. D.

SUPERINTENDENCIA DE INDUSTRIA Y COMERCIO

Rad: 06-031464- -00001-0000 Fea: 2005-04-24 14:08:20

Dep. 2020 NUECREACIONES

Tra: 2 PATENTES Eya 1 REGDEPOSITO

ArJ 328 CTINFORMACION Folios: 44

Ref.: Expediente No. 06-031.464
Solicitud de Patente a nombre de
Ethicon, Inc.

Yo, JIMENA ESCOBAR URIBE, identificada como aparece al pie de mi firma, domiciliada en la Carrera 11 A No. 94-76, of. 305, obrando como apoderada de la sociedad Ethicon, Inc., domiciliada en Somerville, New Jersey, Estados Unidos de América, solicitante de la patente de la referencia, atentamente me permito presentar copia certificada de la solicitud de patente número 11/093,529 presentada en Estados Unidos de América el día 30 de marzo de 2005, con la debida traducción de lo pertinente.

Señor Superintendente,


JIMENA ESCOBAR URIBE
C.C. 39.779.452 de Usaquen.
T.P. 68.228

54/2

ESTADOS UNIDOS DE AMERICA
A TODOS AQUELLOS A QUIENES EL PRESENTE LLEGUE:

DEPARTAMENTO DE COMERCIO DE LOS ESTADOS
UNIDOS

Oficina de Patentes y Marcas de los Estados Unidos

12 de abril de 2006

La presente es para certificar que el anexo es una copia fiel de los registros de la Oficina de Patentes y Marcas de los Estados Unidos, de los documentos de la solicitud de patente identificada más adelante, que cumplieron con los requerimientos para concederse una fecha de presentación de acuerdo con 35 USC 111.

Número de Solicitud: 11/093,529

Fecha de Presentación: 30 de marzo de 2001

Por la autoridad del
Comisionado de Patentes y Marca

EL CÓDIGO DEL PAÍS Y EL NÚMERO DE SU SOLICITUD DE
PRIORIDAD PARA SER USADO EN UNA PRESENTACIÓN
EXTRANJERA BAJO EL CONVENIO DE PARIS, ES US
11/093,529.

(fdo) L. EDELEN
Oficial Certificador
Sello oficial de la Oficina de Patentes y Marcas de los
Estados Unidos de América

RA 7007035



THE UNITED STATES OF AMERICA

TO ALL TO WHOM THESE PRESENTS SHALL COME:

UNITED STATES DEPARTMENT OF COMMERCE

United States Patent and Trademark Office

April 12, 2006

THIS IS TO CERTIFY THAT ANNEXED HERETO IS A TRUE COPY FROM THE RECORDS OF THE UNITED STATES PATENT AND TRADEMARK OFFICE OF THOSE PAPERS OF THE BELOW IDENTIFIED PATENT APPLICATION THAT MET THE REQUIREMENTS TO BE GRANTED A FILING DATE UNDER 35 USC 111.

APPLICATION NUMBER: 11/093,529

FILING DATE: March 30, 2005


THE COUNTRY CODE AND NUMBER OF YOUR PRIORITY APPLICATION, TO BE USED FOR FILING ABROAD UNDER THE PARIS CONVENTION, IS US11/093,529

By Authority of the

**Under Secretary of Commerce for Intellectual Property
and Director of the United States Patent and Trademark Office**

L. Edelen

**L. EDELEN
Certifying Officer**

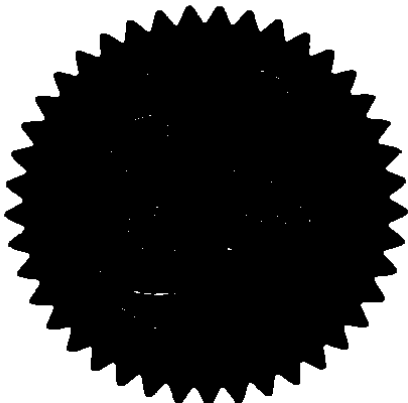

20427 U.S. PTO

UTILITY
PATENT APPLICATION
TRANSMITTAL

Attorney Docket No.

ASP-5025USNP

First Inventor

Ben Fryer, Peter C. Zhu

Title

INTEGRATOR SYSTEM AND METHOD FOR RAPIDLY
DETERMINING EFFECTIVENESS OF A GERMICIDAL
TREATMENT

Express Mail Label No.

EL981642388US

APPLICATION ELEMENTS

See MPEP Chapter 600 concerning utility patent application
contents.

ADDRESS TO:

Commissioner for Patents
P.O. Box 1450
Alexandria, VA 22313-1450

1. ☒ Fee Transmittal Form (e.g., PTO/SB/17)
(submit an original and a duplicate for fee processing)
2. ☐ Applicant claims small entity status.
3. ☒ Specification [Total Pages 26]
(Preferred arrangement set forth below)
- Descriptive Title of the invention
 - Cross Reference to Related Applications
 - Statement Regarding Fed sponsored R&D
 - Reference to sequence listing, a table, or a computer program listing appendix
 - Background of the invention
 - Brief Summary of the invention
 - Brief Description of the Drawings (if filed)
 - Detailed Description
 - Claim(s)
 - Abstract of the Disclosure
4. ☒ Drawing(s) (35 USC 113) [Total Sheets 4]
5. Oath or Declaration [Total Pages 3]
- a. ☒ unexecuted (original or copy)
 - b. ☐ Copy from a prior application (37 CFR 1.63(d))
(for continuation/divisional with Box 18 completed)
 - i. ☐ **DELETION OF INVENTOR(S)**
Signed statement attached deleting
inventor(s) named in the prior application,
see 37 CFR 1.63(d)(2) and 1.33(b).
6. ☐ Application Data Sheet. See 37 CFR 1.76

7. ☐ CD-ROM or CD-R in duplicate, large table or
Computer Program (Appendix)

8. Nucleotide and/or Amino Acid Sequence
Submission (if applicable, all necessary)

- a. ☐ Computer Readable Form (CRF)
- b. ☐ Specification Sequence Listing on:
 - i. ☐ CD-ROM or CD-R (2 copies); or
 - ii. ☐ paper
- c. ☐ Statement verifying identity of above copies

ACCOMPANYING APPLICATION PARTS

9. ☐ Assignment Papers (cover sheet & document(s))
10. ☐ 37 CFR 3.73(b) Statement ☐ Power of Attorney
(when there is an assignee)
11. ☐ English Translation Document (if applicable)
12. ☐ Information Disclosure Statement
(IDS)/PTO-1449 ☐ Copies of IDS Citations
13. ☐ Preliminary Amendment
14. ☒ Return Receipt Postcard (MPEP 503)
(Should be specifically itemized)
15. ☐ Certified Copy of Priority Document(s)
(if foreign priority is claimed)
16. ☐ Request and Certifications under 35 U.S.C. 122
(b)(2)(B)(i). Applicant must attach form
PTO/SB/35 or its equivalent.
17. ☐ Other

18. ☐ If a CONTINUING APPLICATION, check appropriate box and supply the requisite information below and in a
preliminary amendment, or in an Application Data Sheet under 37 CFR 1.76:

☐ Continuation ☐ Divisional ☐ Continuation-in-Part (CIP) of prior application No.: , filed

Prior application information: Examiner Group Art Unit:

For CONTINUATION or DIVISIONAL APPS only: The entire disclosure of the prior application, from which an
oath or declaration is supplied under Box 5b, is considered a part of the disclosure of the accompanying
continuation or divisional application and is hereby incorporated by reference. The incorporation can only be
relied upon when a portion has been inadvertently omitted from the submitted application parts.

19. CORRESPONDENCE ADDRESS

☒ Customer Number or Bar Code Label 000027777 or ☐ Correspondence Address below

Name: Philip S. Johnson, Esq., Johnson & Johnson
One Johnson & Johnson Plaza
New Brunswick, NJ 08933-7003 USA

20. TELEPHONE CONTACT

Please direct all telephone calls or telefaxes to Andrew C. Farmer at:

Telephone: (732) 524-2825 Fax: (732) 524-2808

21. SIGNATURE OF APPLICANT, ATTORNEY, OR AGENT REQUIRED

NAME Andrew C. Farmer Reg. No. 35,868

SIGNATURE

DATE


FEE TRANSMITTAL	<i>Complete If Known</i>	
	Application Number	
	Filing Date	March 30, 2005
	First Named Inventor	Ben Fryer, Peter C. Zhu
	Group Art Unit	
	Attorney Docket Number	ASP-5025USNP

FEE CALCULATION**CLAIMS AS FILED**

(1)	(2)	(3)	(4)	(5)
FOR:	NUMBER FILED	NUMBER EXTRA	RATE	BASIC FEE \$300.00
TOTAL CLAIMS	0 - 20 =	0	x 50.00	\$ 0.00
INDEPENDENT CLAIMS	0 - 3 =	0	x 200.00	\$ 0.00
SEARCH FEE	NA	NA	NA	\$ 500.00
EXAMINATION FEE	NA	NA	NA	\$ 200.00
APPLICATION SIZE FEE If the specification and drawings exceed 100 sheets of paper, the application size fee due is \$250 (\$125 for small entity) for each additional 50 sheets or fraction thereof. See 35 U.S.C. 41(a)(1)(G) and 37 CFR 1.16(a).	0 - 100 =		x 5.00 (for each additional sheet over 100)	\$
MULTIPLE DEPENDENT CLAIMS	☒	N/A	\$360.00	
			TOTAL FEES	\$ 1000.00

METHOD OF PAYMENT

- ☒ Please charge Deposit Account No. 10-0750/ASP-5025USNP/ACF in the amount of \$1,000.00. Three copies of this sheet are enclosed.
- ☒ The Commissioner is hereby authorized to charge any additional fees which may be required in connection with the filing of this communication, or credit any overpayment, to Account No. 10-0750/ASP-5025USNP/ACF. Three copies of this sheet are enclosed.

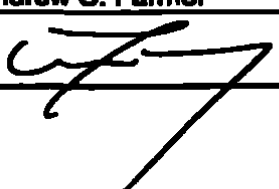SUBMITTED BY:			<i>Complete (if applicable)</i>
Typed or Printed Name	Andrew C. Farmer		Reg. No. 35,868
Signature		Date: March 30, 2005	Deposit Account No. 10-0750

58

DOCKET NO. ASP-5025USNP

IN THE UNITED STATES
PATENT AND TRADEMARK OFFICE

Applicant: Ben Fryer et al.

For : INTEGRATOR SYSTEM AND METHOD FOR RAPIDLY DETERMINING
EFFECTIVENESS OF A GERMIDICAL TREATMENT

Express Mail Certificate

"Express Mail" mailing number: EL961542398US

Date of Deposit: March 30, 2005

I hereby certify that this complete application, including specification pages, claims, drawings and an unexecuted Declaration and Power of Attorney, is being deposited with the United States Postal Service "Express Mail Post Office to Addressee" service under 37 CFR 1.10 on the date indicated above and is addressed to the Commissioner for Patents, P.O. Box 1450, Alexandria, VA 22313-1450.

Laurie Phillips
(Typed or printed name of person mailing paper or fee)

/Laurie Phillips/
(Signature of person mailing paper or fee)

59

INTEGRATOR SYSTEM AND METHOD FOR RAPIDLY DETERMINING EFFECTIVENESS OF A GERMICIDAL TREATMENT

Background of the Invention

Field of the Invention

This invention relates to an integrator system and a method for rapidly determining the effectiveness of a germicidal process for medical equipment.

5 Background

Medical devices are sterilized before use in hospitals, physicians' offices, and other medical facilities. Steam, heat, ethylene oxide, and hydrogen peroxide are commonly used as sterilizing agents.

10 It is standard practice to include a sterility indicator in a load of articles to be sterilized in a sterilizer. The sterility indicator provides a measure of whether the sterilization process was effective in sterilizing the articles in a particular load. If the sterilization process was not effective, as indicated by the sterility indicator, the load of equipment is rejected for use.

15 Biological indicators are generally recognized as reliable sterility indicators. The biological indicator includes a carrier that has been inoculated with spores or other microorganisms. Spores are generally utilized in biological indicators, because spores are more resistant to sterilization than other microorganisms.

20 The biological indicator is placed into the sterilizer with the equipment to be sterilized. At the end of the sterilization process, the biological indicator is removed from the sterilizer, and the carrier is immersed in a sterile culture medium. The culture medium and carrier are incubated for a predetermined time at an appropriate temperature. At the end of the incubation period, it is determined whether any microorganisms have grown in the growth medium. If there is no growth of microorganisms in the growth medium, it is assumed that the equipment in the sterilizer has been properly sterilized. If microorganism
25 growth is observed, the sterilization process was not effective, and the articles in the sterilizer are rejected for use. The growth of microorganisms is determined through a signal such as the generation of turbidity or a color change in a pH indicator due to a pH

60 8

change from byproducts of cell growth in the medium. Biological indicators are described, for example, in U. S. Patent Nos. 5,552,320 and 6,436,659, both of which are incorporated herein by reference in their entirety.

5 Although biological indicators are accurate indicators for the effectiveness of the sterilization cycle, at least 24-48 hours are required to obtain results from the biological indicators. The equipment that was exposed to the sterilization procedure is normally kept in quarantine until the results from the biological indicator are available. Medical equipment is expensive, and storage space in medical facilities is limited. Some hospitals therefore use the equipment before the results are available. Storing quarantined medical
10 equipment is not an efficient use of resources. There is a need for a rapid test for determining the effectiveness of the sterilization process.

Foltz et al. (U. S. Patent No. 6,355,448) describe a method of determining the effectiveness of a sterilization process by using the activity of enzymes rather than spores. It is stated that the enzyme test procedure requires only a few minutes rather than the
15 several days that are required to obtain the results from biological indicators.

The use of a plurality of enzymes rather than a single enzyme was disclosed in U.S. Patent Nos. 5,486,459 and 6,528,277. Using a plurality of enzymes was believed to better mimic the response of a microorganism than a single enzyme.

20 There is a need for sterilization indicators that provide sterilization results rapidly.

SUMMARY OF THE INVENTION

One aspect of the invention involves a method for rapidly determining the effectiveness of an oxidative germicidal process. The method includes providing a substrate having a known amount of a first chemical on the substrate, where the first
25 chemical is selected from the group consisting of a primary amine, mixtures of primary amines, an aldehyde, and mixtures of aldehydes. The first chemical has a first color. The method also includes exposing the substrate and the first chemical to an oxidative germicide, thereby decreasing the known amount of the first chemical to a final amount. The substrate having the final amount of the first chemical having the first color is
30 contacted with a second chemical having a second color, thereby generating a third chemical having a third color. The intensity of the third color is related to the final amount

of the first chemical on the substrate. The second chemical is a chemical selected from the group consisting of a primary amine and mixtures of primary amines when the first chemical is a chemical selected from the group consisting of an aldehyde and a mixture of aldehydes. The second chemical is a chemical selected from the group consisting of an aldehyde and a mixture of aldehydes when the first chemical is a chemical selected from the group consisting of a primary amine and a mixture of primary amines. The method also includes determining the intensity of the third color and determining the effectiveness of the germicidal process from the intensity of the third color.

Advantageously, the effectiveness of the germicidal process is determined by correlating the intensity of the third color with results from biological indicators. In an embodiment, the oxidative germicide is a sterilant. In an alternative embodiment, the oxidative germicide is a disinfectant. The substrate may be an absorbent substrate. Preferably, the substrate is a nonabsorbent substrate.

In an embodiment, the oxidative germicide is a liquid, a vapor, or a gas. Advantageously, the intensity of the third color is determined visually. Preferably, the intensity of said third color is determined spectrophotometrically in the visible or ultraviolet region.

In an embodiment, at least one of the first chemical or the second chemical is colorless. Advantageously, the oxidative germicide is selected from the group consisting of hydrogen peroxide, peracetic acid, ethylene oxide, ozone, and chlorine dioxide. Preferably, the method also includes exposing the substrate and the oxidative germicide to plasma. In an embodiment, the percent completeness of the germicidal process is determined by comparing the intensity of the third color with the intensity of the color of a standard. Preferably, the primary amine is glycine or histidine, and the aldehyde is ortho-phthalaldehyde or glutaldehyde.

Another aspect of the invention involves an integrator for determining the effectiveness of a germicidal process with an oxidative germicide. The integrator includes a substrate with a known amount of a first chemical on the substrate, where the first chemical is selected from the group consisting of a primary amine, mixtures of primary amines, an aldehyde, and mixtures of aldehydes. The substrate is in an enclosure. The first chemical is capable of reacting with the oxidative germicide when exposed to the

62 16

oxidative germicide. The integrator also includes a reservoir of a second chemical, where the second chemical is a chemical selected from the group consisting of a primary amine and mixtures of primary amines when the first chemical is a chemical selected from the group consisting of an aldehyde and a mixture of aldehydes, and the second chemical is a chemical selected from the group consisting of an aldehyde and a mixture of aldehydes when the first chemical is a chemical selected from the group consisting of a primary amine and a mixture of primary amines. The second chemical is capable of reacting with the first chemical to form a third chemical having a color. The reservoir has a breakable barrier that isolates the second chemical from the first chemical and from the oxidative germicide during the contacting of the first chemical with the oxidative germicide. Breaking the breakable barrier in the reservoir contacts the second chemical with the first chemical, thereby forming the third chemical having the color. The reservoir is in the enclosure.

In an embodiment, the breakable barrier in the reservoir includes a frangible ampoule in the enclosure. Advantageously, the integrator also includes a second barrier, where the second barrier is inside the enclosure between the frangible ampoule and the first chemical. The second barrier in the enclosure is permeable to the second chemical. The second barrier prevents fragments from the frangible ampoule from contacting the first chemical.

In an embodiment, the integrator also includes a window in the enclosure, where the window is permeable to the oxidative germicide. The window allows the oxidative germicide to enter the enclosure. Advantageously, the primary amine is selected from the group consisting of glycine and histidine, and the aldehyde is selected from the group consisting of ortho-phthalaldehyde and glutaldehyde.

Preferably, the enclosure on the integrator also includes a transparent window, where the color change on the substrate can be observed through the transparent window visually or with a spectrophotometer.

BRIEF DESCRIPTION OF THE FIGURES

FIG. 1 is a schematic diagram of an integrator according to an embodiment of the present invention;

FIG. 2 is a schematic diagram of a compressible integrator system containing an integrator according to an embodiment of the present invention;

FIG. 3 is a schematic diagram of a slidable integrator system containing an integrator according to an embodiment of the present invention; and

FIG. 4 is a schematic diagram of the slidable integrator system of Figure 3 after an outer shell has been moved over an inner shell of a closable sliding container on the slidable integrator system.

DETAILED DESCRIPTION OF THE INVENTION

The term germicide as used herein is meant to include both sterilants and disinfectants. The term germicidal process as used herein includes both sterilization processes and disinfection processes. Sterilization indicators that utilize chemicals to mimic the resistance of a biological indicator (BI) have been termed as integrators. Integrators utilize an indicator chemical that responds to the germicide that is used in the germicidal process. The chemical reacts with the germicide in a repeatable manner and responds to the factors that are important to sterilization or disinfection in the germicidal process. The reaction of the indicator chemical with the germicide is integrated over time, and the amount of indicator chemical remaining on the integrator is correlated with the BI response.

Integrators integrate the reaction of the chemical over time in response to the critical parameters over a specified range of sterilization cycles.

The integrators and the method according to embodiments of the present invention provide results quickly, reproducibly, and accurately. The chemicals that are used in the integrators are inexpensive and stable. The results that are obtained from the integrators and the method according to embodiments of the present invention correlate well with the results from biological indicator tests.

The integrator according to embodiments of the present invention is meant to mimic the resistance of a biological indicator (BI) without using spores. The integrator

64

according to an embodiment of the present invention includes an indicator chemical that reacts with an oxidative germicide. The integrator is suitable for oxidative germicides including hydrogen peroxide, peracetic acid, ethylene oxide, ozone, and chlorine dioxide. The oxidative germicides may be in the form of a liquid, a vapor, or a gas.

5 In an embodiment, plasma may be utilized in combination with the oxidative germicides to enhance the reaction of the oxidative germicides with the microorganisms in the chamber and the indicator chemical on the integrator and/or to break down the oxidative germicide after use. The use of plasma is optional.

10 The results from the integrator according to embodiments of the present invention are available quickly, approximately 30 seconds to approximately 5-6 minutes, depending on the chemicals selected for the integrators, compared to the 24-48 hours that are normally required to obtain results from a biological indicator.

15 Although described in the context of sterilization with a combination of hydrogen peroxide and plasma with the STERRAD® process, commercially available from Advanced Sterilization Products of Irvine, California, the integrator according to embodiments of the present invention may be used with a variety of germicidal processes. The description of germicidal processes such as sterilization or disinfection with hydrogen peroxide and plasma through the STERRAD® process is illustrative only and is not meant to be limiting.

20 The integrator according to embodiments of the present invention contains an indicator chemical. The indicator chemical reacts with the germicide and responds to the factors that are important for sterilization. The amount of indicator chemical that remains on the integrator after exposure to the germicide can be correlated with the response of BI's that are placed into the sterilization chamber together with the integrator. The
25 response of a BI is a generally accepted measure of the effectiveness of a germicidal process. The response of a BI in the sterilization chamber can be correlated with response of the integrators to "calibrate" the response of the integrators with the response of a biological indicator.

30 In an embodiment, primary amines or aldehydes are used as indicator chemicals in integrators according to embodiments of the present invention. Oxidative germicides react

6513

with both primary amines and aldehydes. Both primary amines and aldehydes are suitable indicator chemicals for integrators according to embodiments of the present invention.

The amount of the primary amine indicator chemical or the aldehyde indicator chemical that remains on the integrator after the germicide process can be used to
5 determine the effectiveness of the sterilization or disinfection process for the load that is treated in the sterilizer.

The amount of primary amine indicator chemical or aldehyde indicator chemical that remains on the integrator can be measured in a variety of ways such as instrumental methods, chemical analysis, etc. Any suitable method of measuring the concentration of
10 the primary amine indicator chemical or the aldehyde indicator chemical is suitable.

In an embodiment, the completeness of the germicidal process may be conveniently determined by observing a change in color in the integrator.

Many primary amines react with aldehydes to form colored products. The amount of primary amine indicator chemical or the amount of aldehyde indicator chemical
15 remaining on the integrator after exposure to the oxidative germicide can be determined from the intensity of the color of the product of the reaction of an aldehyde with a primary amine.

As used herein, an aldehyde that is contacted with a primary amine indicator chemical or a primary amine that is contacted with an aldehyde indicator chemical is
20 termed a "dye precursor", because the product of the reaction of a primary amine with an aldehyde is colored, a "dye", even if neither the primary amine nor the aldehyde has a color.

The primary amine and the aldehyde can change roles, depending on which chemical is utilized as the indicator chemical in the integrator. In an embodiment in which
25 a primary amine is the indicator chemical, the dye precursor is an aldehyde. In an embodiment in which an aldehyde is the indicator chemical, the dye precursor is a primary amine.

The intensity of the color of the colored product resulting from the reaction of the primary amine with the aldehyde can be used to determine the effectiveness of the
30 treatment of the load with the oxidative germicide.

In an embodiment, an integrator according to an embodiment of the present invention contains a first chemical having a first color, where the first chemical is an indicator chemical. The indicator chemical is selected from the group consisting of a primary amine, a mixture of primary amines, an aldehyde, and a mixture of aldehydes.

5 The integrator containing the indicator chemical is placed into a sterilizer with a load of equipment that is to be treated. The load and integrator are contacted with an oxidative germicide in a sterilizer. The oxidative germicide reacts with the indicator chemical, decreasing the amount of indicator chemical remaining on the integrator. The amount of indicator chemical remaining on the integrator after the contacting with the
10 oxidative germicide is a measure of the effectiveness of the germicidal treatment with the oxidative germicide.

At a point when the effectiveness of the treatment with the oxidative germicide is to be determined, the integrator containing the first chemical having the first color is contacted with a second chemical having a second color. The first chemical having the
15 first color is the indicator chemical. The second chemical having the second color is the dye precursor. The dye precursor is a chemical selected from the group consisting of a primary amine, a mixture of primary amines, an aldehyde, and a mixture of aldehydes. Primary amines may not be mixed with aldehydes to form the dye precursor. In an
embodiment where the first chemical, the indicator chemical, is a primary amine or a
20 mixture of primary amines, the second chemical, the dye precursor, is an aldehyde or a mixture of aldehydes. In an embodiment where the first chemical, the indicator chemical, is an aldehyde or a mixture of aldehydes, the second chemical, the dye precursor, is a primary amine or a mixture of primary amines.

The product from the reaction of the first chemical, the indicator chemical, with the
25 second chemical, the dye precursor, is a third chemical having a third color. The intensity of the third color of the third chemical resulting from the reaction of the first chemical with the second chemical can be used to determine how much of the first chemical, the indicator chemical, remains on the integrator. The amount of the first chemical indicator chemical that remains on the integrator is a measure of how effective the treatment with the
30 oxidative germicide was. If only a small amount of the first chemical indicator chemical remains on the integrator, the intensity of the third color due to the third chemical is low.

67 15

A low intensity of the third color is an indication that the treatment with the oxidative germicide was effective.

In an embodiment, the degree of completeness of the treatment with the oxidative germicide can be determined from the intensity of the third color due to the third chemical, the product of the reaction of the first chemical indicator compound with the second chemical dye precursor. The indicator compound on the integrator, the first compound, reacts with the oxidative germicide in parallel with the germicidal treatment of the load in the sterilization chamber. The intensity of the third color due to the third chemical resulting from the reaction of the first chemical indicator chemical with the second chemical dye precursor decreases as the amount of the indicator chemical decreases due to the reaction with the oxidative germicide.

The intensity of the third color can be correlated with results from biological indicators that are placed into the sterilization chamber together with the indicators. The intensity of the third color can be correlated with the percent sterilization or percent disinfection, as determined from biological indicators. The percent sterilization or disinfection can therefore be determined from the intensity of the third color due to the third compound.

In an embodiment, the first chemical, the indicator chemical, is placed on a substrate for ease of handling. The substrate can be a variety of materials. The substrate can be an absorbing substrate or a nonabsorbing substrate. Absorbing substrates absorb the germicide. Nonabsorbing substrates absorb little or none of the germicide.

Filter paper is an absorbing substrate, because filter paper absorbs the germicide. A glass filter disk is a nonabsorbing substrate, because the glass filter disk does not absorb significant quantities of the germicide and is thus preferred.

The indicator chemical can be packaged in a water-soluble binder such as an acrylic polymer or carboxymethylcellulose. The indicator chemical and water-soluble binder can be applied to the surface of the integrator or sterility indicator by, for example, ink jet printing a solution of the indicator chemical and the water-soluble binder onto the surface of an inert backing material.

Absorbing substrates absorb germicide during the germicidal process. When the second chemical dye precursor is contacted with the absorbing substrate, the absorbed

68 12

germicide on the absorbing substrate can react with the dye precursor. Oxidative germicides generally react with primary amines and aldehydes, the two forms of dye precursor. It is therefore generally advantageous to use excess second chemical primary amine or aldehyde dye precursor when the substrate is an absorbing substrate, because the absorbed germicide reacts with the dye precursor when the dye precursor is contacted with the absorbing substrate.

In an embodiment, the integrator containing the indicator chemical is placed in the sterilizer with the equipment to be sterilized and is exposed to the germicide. The indicator chemical reacts with the germicide, reducing the initial concentration of the indicator chemical from an initial value to a final value.

After the germicidal process is complete, the integrator containing the first chemical indicator chemical is exposed to the second chemical dye precursor. If any indicator chemical is still present, contacting the indicator chemical with the dye precursor forms the third compound having the third color. If significant color develops on the integrator, the germicidal cycle is judged to have not been effective.

It is generally preferred that the second chemical dye precursor be contacted with the integrator after the conclusion of the cycle, because the dye precursor reacts with the germicide. If the dye precursor is contacted with the integrator before the conclusion of the cycle, the germicide will react with the second chemical dye precursor, and it may be necessary to add dye precursor to provide sufficient dye precursor to cause a color change from the reaction of the second chemical dye precursor with the first chemical indicator chemical, forming the third chemical having the third color. The second chemical dye precursor is therefore generally contacted with the integrator at the end of the cycle. In an embodiment, the cycle may be a canceled cycle.

In an embodiment, the second chemical is isolated from the germicide until the conclusion of the cycle. Isolating the second chemical dye precursor from the germicide protects the second chemical from reacting with the germicide and being destroyed.


The color change from the reaction of the first chemical indicator chemical with the second chemical dye precursor to form the third chemical having the third color can be determined visually. Because a visual change is somewhat subjective, the color change is generally determined with an optical detector. The optical detector for the color change

resulting from the reaction of the first chemical indicator chemical with the second chemical dye precursor can operate at visible or ultraviolet wavelengths.

The primary amine may be any suitable primary amine. In an embodiment, the primary amine is an amino acid. In an embodiment, the primary amine is selected from the group consisting of arginine, histidine, and combinations thereof. Other suitable primary amines include the following amino acids: alanine, proline, amino-caproic acid, phenylalanine, tryptophan, methionine, glycine, serine, cysteine, tyrosine, glutamine, aspartic acid, glutamic acid, lysine, arginine, and histidine. Peptides or polypeptides formed from any number or type of amino acids are also suitable primary amines.

Arginine is an exemplary primary amine indicator chemical. Arginine gives a strong, rapid color change when exposed to aldehydes. Arginine also reacts rapidly with germicides. Arginine is a water-soluble solid that is conveniently weighed, dissolved in a solvent, and applied to the substrate or other support. Other primary amines can be used in other embodiments, as will become clear in the description and the Examples below.

Arginine has structure I, below.

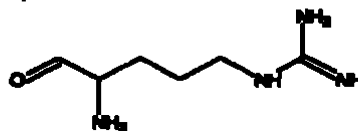

I

The NH_2 groups are primary amine groups. The NH groups are secondary amine groups. Aldehydes often do not react with secondary amine groups.

The aldehyde may be any aldehyde that reacts with primary amines but not secondary or tertiary amines to produce a color. Aldehydes such as OPA (ortho-phthalaldehyde), glutaldehyde, and aromatic aldehydes are suitable. Other aldehydes are also suitable.

Figure 1 shows a diagram of an integrator system 10 according to an embodiment of the present invention. The integrator system 10 of Figure 1 includes integrator chemistry 14 located on integrator strip 16, where the integrator strip 16 is an inert material for supporting the integrator chemistry 14. The integrator strip 16 is generally made of a

material that does not react with or adsorb the germicide. The integrator chemistry 14 includes a first chemical, the indicator chemical. The integrator strip 16 is a substrate for the integrator chemistry 14.

Peel off label 18 is optionally located on the integrator strip 16. Information on the sterilization cycle can be written on the peel off label 18, and the peel off label 18 with the information on the sterilization cycle can be placed into a sterilization logbook. Chemical indicator strip 20 contains a chemical that undergoes a color change when exposed to the germicide. A color change in chemical indicator strip 20 simply shows that the chemical indicator strip 20 has been exposed to the germicide. The chemical indicator strip 20 is not an indicator of the effectiveness of sterilization but is simply an indicator as to whether the chemical indicator strip 20 has been exposed to germicide.

The color change on the chemical indicator strip 20 shows the operator that the integrator system 10 should not be used again. Bordeaux Red changes color when exposed to hydrogen peroxide. Other dyes can be used on the chemical indicator strip 20 to indicate exposure to other germicides. Suitable dyes are described, for example, in U. S. Patent No. 5,942,438, which is incorporated herein by reference in its entirety.

After the cycle, the integrator chemistry 14 portion of the integrator strip 16 is exposed to the second chemical, the dye precursor. The second chemical dye precursor reacts with the first chemical indicator chemical on the integrator chemistry 14 to form the third chemical having the third color. The presence of a significant amount of color on the integrator chemistry 14 on integrator strip 16 after exposure of the integrator chemistry 14 to the second chemical dye precursor indicates that the cycle was not effective.

Figure 2 shows a compressible integrator system 22. Compressible integrator system 22 of Figure 2 includes integrator chemistry 14 on substrate 44 located in container 26. Gas permeable surface 24 allows the germicide to enter container 26 and contact the integrator chemistry 14.

The dye precursor 28 is contained in reservoir 30. Supports 32 are located adjacent to reservoir 30. Reservoir 30 protects the dye precursor from being destroyed by reacting with the oxidative germicide during the germicidal cycle.

After the cycle, compressible integrator system 22 is crushed or squeezed. Supports 32 pierce the reservoir 30, and the second chemical dye precursor 28 that was

contained in the reservoir 30 contacts the integrator chemistry 14. The second chemical dye precursor reacts with any first chemical indicator chemical that remains after the cycle. If any first chemical indicator chemical remains on the substrate 44, the first chemical indicator chemical in the integrator chemistry 14 reacts with the second chemical dye precursor to form the third chemical having the third color on the substrate 44, indicating that the germicidal process was incomplete. A lack of color on the integrator chemistry 14 indicates that the germicidal treatment was successful.

Figure 3 shows a schematic diagram of a slidable integrator system 34. Slidable integrator system 34 includes closable sliding container 36. The closable sliding container 36 is formed of outer shell 38 and inner shell 40. Outer shell 38 slides over inner shell 40. Sliding outer shell 38 over inner shell 40 opens window 42 in closable sliding container 36. Window 42 allows germicide to enter the interior of closable sliding container 36. The germicide may be a liquid, a vapor, or a gas.

The closable sliding container 36 contains substrate 44 supporting integrator chemistry 14. The integrator chemistry 14 includes the first chemical, the indicator chemical. Substrate 44 is located adjacent transparent window 46 on inner shell 40. Any color change in the substrate 44 can be observed through transparent window 46. Substrate 44 is the substrate for the indicator chemistry 14. The indicator chemistry 14 includes the first chemical, the indicator chemical. In an embodiment, substrate 44 is a glass filter, a nonabsorbing substrate.

The second chemical, the dye precursor 24, is contained in crushable ampoule 48 inside closable sliding container 36. Crushable ampoule 48 is made of a frangible material such as glass. Crushable ampoule 48 protects the second chemical, the dye precursor 24, from being destroyed by the germicide during the germicidal process.

Wedge 50 is attached to an inside of outer shell 38 of the closable sliding container 36. Wedge 50 is a projection on the inside of the outer shell 38. In an embodiment wedge 50 has a sharp edge to aid in penetrating the crushable ampoule 48. Barrier 52 is located inside closable sliding container 34 between crushable ampoule 48 and substrate 44. Barrier 52 prevents fragments of crushable ampoule 48 from interfering with the reading of substrate 44. Barrier 52 is permeable to the second chemical, the dye precursor 24. When crushable ampoule 48 is crushed, the second chemical, the dye precursor, is released and

can flow through barrier 52 to contact the first chemical indicator chemical on substrate 44. In an embodiment, barrier 52 is a wire screen.

Figure 4 shows a schematic diagram of the slidable integrator system 34 and closable sliding container 36 of Figure 3 after the conclusion of the cycle. The outer shell 38 of the closable sliding container 36 in Figure 4 has been moved toward the left of Figure 4 by sliding outer shell 38 over inner shell 40.

Sliding outer shell 38 over inner shell 40 has several effects, as shown in Figure 4. First, sliding outer shell 38 over inner shell 40 closes window 42. Closing window 42 isolates closable sliding container 36, retaining the second chemical dye precursor inside the closable sliding container 36. The second chemical dye precursor can stain the hands of the operator. Second, sliding the outer shell 38 over inner shell 40 forces wedge 50 into contact with crushable ampoule 48, pushing the crushable ampoule 48 into contact with barrier 52, crushing the crushable ampoule 48. Crushing the crushable ampoule 48 releases the second chemical, the dye precursor 24. The second chemical dye precursor 24 that is released flows through barrier 52 and contacts the substrate 44.

If any first chemical indicator chemical remains on the substrate 44 when the second chemical dye precursor 24 contacts the substrate 44, the second chemical dye precursor 24 reacts with the first chemical indicator chemical in the integrator chemistry 14 to form the third chemical having a third color. The third color is distinctive and readily distinguished from the first color of the first chemical and the second color of the second chemical. Barrier 52 prevents fragments of the crushable ampoule 48 from contacting the substrate 44 and interfering with the determination. The color change on the substrate 44 can be observed through transparent window 46.

Method of Using the Integrator

The integrator according to an embodiment of the present invention is placed into the sterilization chamber with the load that is to be sterilized, and the germicidal cycle is run. After the completion of the sterilization cycle, the second chemical dye precursor is contacted with the integrator. The second chemical, the dye precursor, reacts with any first chemical indicator chemical remaining on the integrator to produce the third chemical having the third color. The color that is produced depends on the structure of the first chemical indicator chemical and the second chemical dye precursor. The intensity of the

color depends on the amount of first chemical indicator chemical that remains on the integrator and on the concentration of the second chemical dye precursor.

If color develops on the integrator within a predetermined time period, such as approximately 5-6 minutes after the second chemical dye precursor is added to the integrator, the sterilization cycle is judged to have been ineffective.

The intensity of the color on the integrator may be judged visually by comparing to a color standard, or the intensity of the color may be measured spectrophotometrically in the visible or ultraviolet region. Judging the color intensity visually is more subjective than measuring the color intensity with an instrument.

The second chemical dye precursor may be contacted with the integrator in various manners. In an embodiment, the second chemical dye precursor is contacted with the integrator manually using a pipette, an eyedropper, or other suitable device.

Manual addition of the second chemical dye precursor is appropriate with an integrator such as the integrator shown in Figure 1, where there is no method of protecting the second chemical dye precursor from being exposed to and being destroyed by the germicide during the sterilization cycle.

Figures 2 and 3 show embodiments of integrators where the second chemical dye precursor is present in the integrator during the sterilization cycle but is protected from exposure to the germicide by being enclosed in a reservoir or a crushable ampoule. Other means of protecting the second chemical dye precursor from being exposed to the germicide may be used in other embodiments.

The following examples are meant to be illustrative only and are not meant to be limiting on the scope.

EXAMPLES

Example 1

**Results of Tests with an Integrator Containing Arginine as the Indicator
Chemical With Varying Injection Volumes of Hydrogen Peroxide**

5 A series of integrators was prepared by contacting paper disks with an aqueous solution of arginine. The integrators were placed in a STERRAD® 50 sterilizer with a load of equipment to be sterilized and several biological indicators. The paper disks are absorbing substrates.

10 The sterilizer was evacuated to 0.8 torr. Plasma was produced in the chamber for 15 minutes to condition the load. The sterilizer was evacuated further to 0.4 torr, and hydrogen peroxide was injected and contacted with the load, integrators, and biological indicators for 6 minutes.

15 The sterilizer was vented with air for 2 minutes. The sterilizer was evacuated again to 0.5 torr, and plasma was produced for an additional 2 minutes. The plasma power was 400 watts for both plasma exposures.

The paper integrator disks were contacted with 100 µL of a 5% solution of ortho-phthalaldehyde (OPA) dye precursor after the germicide cycle, and the response of the integrators was measured visually after the integrators had been contacted with the OPA.
20 The results are shown in Table 1 below.

25
75

Table 1
Integrator and Biological Indicator Test Results

Hydrogen Peroxide Injection Volume μL	Biological Indicator Results	Integrator Color Intensity
100	100% positive	Dark orange color in 3 minutes
250	40% positive	Orange color after 5-6 minutes, a few are colorless
400	0% positive	10-15% show faint color after 5-6 minutes
500	0% positive	No color

5 The cycles with 100 and 250 μL of hydrogen peroxide were ineffective, as shown by the positive biological indicator results. The integrator results were consistent with the biological indicator results, because the integrators for the cycles with 100 and 250 μL of hydrogen peroxide had significant color within 3-6 minutes after being contacted with the OPA. The color is the result of a reaction between unreacted arginine, the primary amine indicator compound, and OPA, the aldehyde dye precursor.

10 In the sterilization cycle with 400 μL of hydrogen peroxide, 10-15% of the integrators developed a faint color after 5-6 minutes. None of the biological indicators in this run were positive. The faint color that developed on the integrators after 5-6 minutes is a showing that the integrators provide a more stringent measure of the degree of sterilization than the biological indicators.

15 None of the biological indicators in the run with 500 μL of hydrogen peroxide were positive. None of the integrators had any color. Both the biological indicator results and the integrator results are consistent in showing that the sterilization cycle with 500 μL of hydrogen peroxide was effective.

20 The data from the integrators according to an embodiment of the present invention were consistent with the data from the biological indicators. However, the results from the

24
76

integrators were available in 5-6 minutes, compared to 24-48 hours for the biological indicator results.

Example 2

Integrator Response Measured With a Spectrometer

5 Non adsorbent glass fiber disks were impregnated with an aqueous solution of arginine. The glass fiber disks are nonabsorbing substrates. The disks were placed in a STERRAD® 50 sterilizer and processed under the same conditions as in Example 1 with varying injection volumes of hydrogen peroxide. The quantities of hydrogen peroxide are shown in Table 2 below.

10 The glass fiber disks were contacted with 50 µL of a 5% aqueous solution of OPA dye precursor after the cycle was complete. The intensity of light absorption of the disks was determined with a TAOS TCS230EVM evaluation module color sensor (Parallax, Rocklin, California) at approximately 470 nm, approximately 550 nm, and approximately 610 nm (red, green, and blue wavelengths).

15 Hydrogen peroxide injection volumes of 50, 300, and 1000 µL were used. All of the BI's would be negative with a hydrogen peroxide injection of 300 µL in the STERRAD® 50 sterilizer. Table 2 summarizes the color intensities at the 610 nm (blue) wavelength. The integrator responses shown in Table 2 are the differences between the initial reading and the final reading at 30 seconds after addition of OPA. Large numbers
20 for the integrator response indicate more color and less effective sterilization.

Table 2

Hydrogen Peroxide Injection Volume (µL)	Integrator Response (Range) for Samples
1000	5-17
300	28-47
50	58-85

Small numbers in the integrator response indicate effective sterilization. The range of 5-17 in the integrator response for the samples that were exposed to 1000 µL of
25 hydrogen peroxide is consistent with effective sterilization. The integrator response range

of 28-47 for the samples that were exposed to 300 μ L of hydrogen peroxide is consistent with effective sterilization.

The integrator response range of 58-85 for the samples that were exposed to 50 μ L of hydrogen peroxide indicates significant color, showing ineffective sterilization.

5 The results of Example 2 demonstrate that the results from the integrators according to embodiments of the present invention can be measured effectively with a spectrophotometer rather than visually, as in Example 1.

Example 3

Integrator Tests With Histidine as the Indicator Chemical and OPA as the Dye Precursor

10 Histidine rather than arginine was used as the primary amine indicator chemical to form a series of integrators in Example 3. An aqueous solution of histidine was placed on a series of glass fiber disks to form integrators according to an alternative embodiment of the present invention.

15 The integrators were placed into a standard STERRAD[®] 50 validation load, and cycles were run as described in Examples 1 and 2.

Cycles were run with 0 μ L of hydrogen peroxide and 300 μ L of hydrogen peroxide. The integrators were contacted with 50 μ L of a 5 volume % solution of OPA at the end of the sterilization process. The intensity of the third color of the third compound on the
20 integrators was measured with the TAOS color sensor described in Example 2. A combination of red, green, and blue (RGB) wavelengths (470, 550, and 610 nm) were used to measure the response, because the color of the product of histidine and OPA is different than the color of the product of arginine and OPA. A root mean square (RMS) of the absorption was calculated for all three wavelengths measured by the sensor.

25 The light absorption results were consistent with the treatment with 0 μ L of hydrogen peroxide being ineffective at sterilization and the treatment with 300 μ L of hydrogen peroxide being effective at sterilization. The example demonstrates that histidine can be used as a primary amine indicator compound in place of arginine. A wide variety of primary amines can be used as primary amine indicator compounds.

Example 4

Integrator Tests With Histidine and Glutaraldehyde

5 A series of integrators was prepared with histidine as the indicator chemical on glass fiber disks as the substrate. The integrators were processed in a STERRAD® 50 sterilizer under the same conditions as in Example 3. Glutaraldehyde rather than OPA was contacted with the processed integrators as the aldehyde dye precursor. The results are shown in Table 4 below.

Table 4

Results From Integrator Tests With Histidine and Glutaraldehyde

Volume of Hydrogen Peroxide (μL)	Range in the Change in Red Wavelength Reading (Change = Initial Color Reading – Color Reading at 30 Seconds)
0	11-29
300	0-7

10

The small change in light absorption from the start time to the finish time for the integrators that were exposed to 300 μL of hydrogen peroxide demonstrates that the sterilization with 300 μL of hydrogen peroxide was effective.

15 The large change in light absorption from the start time to the reading at 30 seconds for the integrators that were exposed to a volume of 0 μL of hydrogen peroxide is a showing that the sterilization was not effective.

The results of Example 4 are a demonstration that glutaraldehyde can be used as an aldehyde dye precursor in place of OPA

Example 5

20

Integrators With OPA as the Indicator Chemical

A series of integrators are prepared with OPA as the indicator chemical. The indicators are placed into a sterilizer with a load to be sterilized and a series of biological indicators. The load, biological indicators, and integrators are contacted with 100-500 μL of hydrogen peroxide under the conditions described in Example 1.

28
79

The integrators with OPA as the indicator chemical are contacted with an aqueous solution of arginine as the dye precursor. The results from the integrators with OPA as an aldehyde indicator chemical and arginine as a primary amine dye precursor correlate well with the results from the biological indicators.

5 The results from Example 5 demonstrate that aldehydes such as OPA can be used as the indicator chemical with primary amines such as arginine as the dye precursor.

The integrator according to embodiments of the present invention allows the effectiveness of the sterilization process to be determined rapidly.

10 Various modifications and alterations of this invention will be apparent to those skilled in the art without departing from the scope and spirit of this invention. It is to be understood that the invention is not limited to the embodiments disclosed therein, and that the claims should be interpreted as broadly as the prior art allows.

NOVEDAD DE LA INVENCION

REIVINDICACIONES

- 5 1.- Un método para determinar rápidamente la eficacia de un proceso germicida oxidativo, que comprende: proveer un sustrato que tiene una cantidad conocida de un primer reactivo sobre el sustrato, en donde dicho primer reactivo se selecciona del grupo que consiste de una amina primaria, mezclas de aminas primarias, un aldehído, y mezclas de aldehídos;
- 10 en donde dicho primer reactivo tiene un primer color; exponer el sustrato que tiene la cantidad conocida del primer reactivo a un germicida oxidativo, reduciendo con ello la cantidad conocida del primer reactivo a una cantidad final del primer reactivo; poner en contacto el sustrato que tiene la cantidad final del primer reactivo que tiene el primer color, con un segundo reactivo que
- 15 tiene un segundo color, generando así un tercer reactivo que tiene un tercer color, dicho tercer color teniendo una intensidad, en donde la intensidad de dicho tercer color se relaciona con la cantidad final de dicho primer reactivo sobre dicho sustrato, en donde dicho segundo reactivo es un reactivo seleccionado del grupo que consiste de una amina primaria y mezclas de
- 20 aminas primarias cuando dicho primer reactivo es un reactivo seleccionado del grupo que consiste de un aldehído y una mezcla de aldehídos, y dicho segundo reactivo es un reactivo seleccionado del grupo que consiste de un aldehído y una mezcla de aldehídos cuando dicho primer reactivo es un

21
81

reactivo seleccionado del grupo que consiste de una amina primaria y una mezcla de aminas primarias; determinar la intensidad de dicho tercer color, y determinar la eficacia del proceso germicida a partir de la intensidad de dicho tercer color.

5 2.- El método de conformidad con la reivindicación 1, caracterizado además porque el paso de determinar la eficacia del proceso germicida comprende correlacionar la intensidad de dicho tercer color con resultados de indicadores biológicos.

10 3.- El método de conformidad con la reivindicación 1, caracterizado además porque dicho germicida oxidativo es un esterilizante.

4.- El método de conformidad con la reivindicación 1, caracterizado además porque dicho germicida oxidativo es un desinfectante.

5.- El método de conformidad con la reivindicación 1, caracterizado además porque dicho sustrato es un sustrato absorbente.

15 6.- El método de conformidad con la reivindicación 1, caracterizado además porque dicho sustrato es un sustrato no absorbente.

7.- El método de conformidad con la reivindicación 1, caracterizado además porque dicho germicida oxidativo es un líquido, un vapor, o un gas.

20 8.- El método de conformidad con la reivindicación 1, caracterizado además porque la intensidad de dicho tercer color se determina visualmente.

9.- El método de conformidad con la reivindicación 1,

caracterizado además porque la intensidad de dicho tercer color se determina espectrofotométricamente en la región visible o ultravioleta.

10.- El método de conformidad con la reivindicación 1, caracterizado además porque por lo menos uno de dicho primer reactivo y
5 dicho segundo reactivo, es incoloro.

11.- El método de conformidad con la reivindicación 1, caracterizado además porque dicho germicida oxidativo se selecciona del grupo que consiste de peróxido de hidrógeno, ácido peracético, óxido de etileno, ozono y dióxido de cloro.

10 12.- El método de conformidad con la reivindicación 1, caracterizado además porque comprende exponer dicho sustrato y dicho germicida oxidativo a plasma.

13.- El método de conformidad con la reivindicación 1, caracterizado además porque se determina un porcentaje de consumación del
15 proceso germicida comparando la intensidad del tercer color con la intensidad de color de un estándar.

14.- El método de conformidad con la reivindicación 1, caracterizado además porque la amina primaria se selecciona del grupo que consiste de glicina e histidina, y el aldehído se selecciona del grupo que
20 consiste de orto-ftalaldehído y glutaldehído.

15.- Un integrador para determinar la eficacia de un proceso germicida con un germicida oxidativo, dicho integrador comprendiendo: un sustrato con una cantidad conocida de un primer reactivo sobre el sustrato,

en donde dicho primer reactivo se selecciona del grupo que consiste de una amina primaria, mezclas de aminas primarias, un aldehído, y mezclas de aldehídos; en donde dicho sustrato está en un recinto; en donde dicho primer reactivo es capaz de reaccionar con dicho germicida oxidativo cuando se expone a dicho germicida oxidativo; un depósito de un segundo reactivo, en donde dicho segundo reactivo es un reactivo seleccionado del grupo que consiste de una amina primaria y mezclas de aminas primarias cuando dicho primer reactivo es un reactivo seleccionado del grupo que consiste de un aldehído y una mezcla de aldehídos, y dicho segundo reactivo es un reactivo seleccionado del grupo que consiste de un aldehído y una mezcla de aldehídos cuando dicho primer reactivo es un reactivo seleccionado del grupo que consiste de una amina primaria y una mezcla de aminas primarias; en donde dicho segundo reactivo es capaz de reaccionar con dicho primer reactivo para formar un tercer reactivo que tiene un color; en donde dicho depósito tiene una barrera rompible que aísla el segundo reactivo del primer reactivo y del germicida oxidativo durante el contacto del primer reactivo con el germicida oxidativo, en donde el rompimiento de dicha barrera rompible del depósito pone en contacto el segundo reactivo con el primer reactivo, formando así dicho tercer reactivo que tiene el color; y en donde el depósito está en dicho recinto.

16.- El integrador de conformidad con la reivindicación 15, caracterizado además porque dicha barrera rompible de dicho depósito comprende una ampolla frangible en dicho recinto.

37
84

17.- El integrador de conformidad con la reivindicación 16, caracterizado además porque comprende una segunda barrera, en donde dicha segunda barrera está dentro de dicho recinto, entre dicha ampolla frangible y dicho primer reactivo, en donde dicha segunda barrera de dicho
5 recinto es permeable al segundo reactivo, y en donde dicha segunda barrera impide que fragmentos de la ampolla frangible hagan contacto con el primer reactivo.

18.- El integrador de conformidad con la reivindicación 15, caracterizado además porque comprende una ventana en dicho recinto, en
10 donde dicha ventana es permeable a dicho germicida oxidativo, dicha ventana permitiendo que dicho germicida oxidativo entre a dicho recinto.

19.- El integrador de conformidad con la reivindicación 15, caracterizado además porque la amina primaria se selecciona del grupo que consiste de glicina e histidina, y el aldehído se selecciona del grupo que
15 consiste de orto-ftalaldehído y glutaldehído.

20.- El integrador de conformidad con la reivindicación 15, caracterizado además porque dicho recinto también comprende una ventana transparente, en donde un cambio de color en dicho sustrato puede ser observado a través de dicha ventana transparente, visualmente o con un
20 espectrofotómetro.

WHAT IS CLAIMED IS:

1. A method for rapidly determining the effectiveness of an oxidative germicidal process, said method comprising:
providing a substrate having a known amount of a first chemical on the substrate, wherein said first chemical is selected from the group consisting of a primary amine, mixtures of primary amines, an aldehyde, and mixtures of aldehydes, wherein said first chemical has a first color;
exposing the substrate having the known amount of the first chemical to an oxidative germicide, thereby decreasing the known amount of the first chemical to a final amount of the first chemical;
contacting the substrate having the final amount of the first chemical having the first color with a second chemical having a second color, thereby generating a third chemical having a third color, said third color having an intensity, wherein the intensity of said third color is related to the final amount of said first chemical on said substrate, wherein said second chemical is a chemical selected from the group consisting of a primary amine and mixtures of primary amines when said first chemical is a chemical selected from the group consisting of an aldehyde and a mixture of aldehydes, and said second chemical is a chemical selected from the group consisting of an aldehyde and a mixture of aldehydes when said first chemical is a chemical selected from the group consisting of a primary amine and a mixture of primary amines, ;
determining the intensity of said third color; and
determining the effectiveness of the germicidal process from the intensity of said third color.
2. The method of claim 1, wherein the step of determining the effectiveness of the germicidal process comprises correlating the intensity of said third color with results from biological indicators.
3. The method of claim 1, wherein said oxidative germicide is a sterilant.

7A
86

4. The method of claim 1, wherein said oxidative germicide is a disinfectant.
5. The method of claim 1, wherein said substrate is an absorbent substrate.
- 5 6. The method of claim 1, wherein said substrate is a nonabsorbent substrate.
7. The method of claim 1, wherein said oxidative germicide is a liquid, a vapor, or a gas.
- 10 8. The method of claim 1, wherein the intensity of said third color is determined visually.
9. The method of claim 1, wherein the intensity of said third color is determined spectrophotometrically in the visible or ultraviolet region.
- 15 10. The method of claim 1, wherein at least one of said first chemical and said second chemical is colorless.
11. The method of claim 1, wherein said oxidative germicide is selected from the group consisting of hydrogen peroxide, peracetic acid, ethylene oxide, ozone, and chlorine dioxide.
- 20 12. The method of claim 1, further comprising exposing said substrate and said oxidative germicide to plasma.
- 25 13. The method of claim 1, wherein a percent completeness of the germicidal process is determined by comparing the intensity of the third color with an intensity of a color of a standard.

30

38
67

14. The method of claim 1 wherein the primary amine is selected from the group consisting of glycine and histidine and the aldehyde is selected from the group consisting of ortho-phthalaldehyde and glutaldehyde.

5 15. An integrator for determining the effectiveness of a germicidal process with an oxidative germicide, said integrator comprising:

a substrate with a known amount of a first chemical on the substrate, wherein said first chemical is selected from the group consisting of a primary amine, mixtures of primary amines, an aldehyde, and mixtures of aldehydes;

10 wherein said substrate is in a enclosure;

wherein said first chemical is capable of reacting with said oxidative germicide when exposed to said oxidative germicide;

15 a reservoir of a second chemical, wherein said second chemical is a chemical selected from the group consisting of a primary amine and mixtures of primary amines when said first chemical is a chemical selected from the group consisting of an aldehyde and a mixture of aldehydes, and said second chemical is a chemical selected from the group consisting of an aldehyde and a mixture of aldehydes when said first chemical is a chemical selected from the group consisting of a primary amine and a mixture of primary amines;

20 wherein said second chemical is capable of reacting with said first chemical to form a third chemical having a color;

25 wherein said reservoir has a breakable barrier that isolates the second chemical from the first chemical and from the oxidative germicide during the contacting of the first chemical with the oxidative germicide, wherein breaking said breakable barrier in the reservoir contacts said second chemical with said first chemical, thereby forming said third chemical having the color; and

30 wherein said reservoir is in said enclosure.

36
48

- 5
- 10
- 15
- 20
- 16. The integrator of claim 15, wherein said breakable barrier in said reservoir comprises a frangible ampoule in said enclosure.
 - 17. The integrator of claim 16, further comprising a second barrier, wherein said second barrier is inside said enclosure between said frangible ampoule and said first chemical, wherein said second barrier in said enclosure is permeable to said second chemical, and wherein said second barrier prevents fragments from the frangible ampoule from contacting the first chemical.
 - 18. The integrator of claim 15, further comprising a window in said enclosure, wherein said window is permeable to said oxidative germicide, said window allowing said oxidative germicide to enter said enclosure.
 - 19. The integrator of claim 15, wherein the primary amine is selected from the group consisting of glycine and histidine and the aldehyde is selected from the group consisting of ortho-phthalaldehyde and glutaldehyde.
 - 20. The integrator of claim 15, wherein said enclosure further comprises a transparent window, wherein a color change on said substrate can be observed through said transparent window visually or with a spectrophotometer.

Abstract

**INTEGRATOR SYSTEM AND METHOD FOR RAPIDLY
DETERMINING EFFECTIVENESS OF A GERMICIDAL TREATMENT**

- 5 The effectiveness of an oxidative sterilization process is determined by exposing a substrate having a known amount of a primary amine or aldehyde indicator chemical to an oxidative germicide. The oxidative germicide reacts with the indicator chemical. The amount of indicator chemical remaining after exposure to the germicide is determined by reacting the indicator chemical with a dye precursor chemical to form a colored product.
- 10 The amount of indicator chemical remaining on the substrate is determined from the intensity of the color of the colored product. The amount of indicator chemical remaining on the substrate is a measure of the effectiveness of the germicidal treatment. The dye precursor is an aldehyde when the indicator chemical is a primary amine and a primary amine when the indicator chemical is an aldehyde. An integrator for determining the
- 15 effectiveness of the germicidal process includes a substrate and an indicator chemical, where the indicator chemical is a primary amine or an aldehyde.

Please type a plus sign (+) inside this box ☐

PTO/SB/01 (10-03)

Approved for use through 10/31/2002. OMB 0851-0032

U.S. Patent and Trademark Office; U.S. DEPARTMENT OF COMMERCE

Under the Paperwork Reduction Act of 1995, no person is required to respond to a collection of information unless it contains a valid OMB control number.

DECLARATION AND POWER OF ATTORNEY FOR UTILITY OR DESIGN PATENT APPLICATION (37 CFR 1.63)		Attorney Docket Number	ASP-5025USNP	
☒ Declaration Submitted with Initial Filing OR ☐ Declaration Submitted after Initial Filing (Surcharge (37 CFR 1.16(e)) required)		First Named Inventor	Ben Fryer et al.	
		COMPLETE IF KNOWN		
		Application Number		
		Filing Date	March 30, 2005	
		Group Art Unit		
		Examiner Name		
As a below named inventor, I hereby declare that: My residence, mailing address, and citizenship are as stated below next to my name. I believe I am the original, first and sole inventor (if only one name is listed below) or an original, first and joint inventor (if plural names are listed below) of the subject matter which is claimed and for which a patent is sought on the invention entitled: INTEGRATOR SYSTEM AND METHOD FOR RAPIDLY DETERMINING EFFECTIVENESS OF A GERMICIDAL TREATMENT (Title of the Invention) the specification of which ☒ is attached hereto OR ☐ was filed on as United States Application Number or PCT International Application Number and was amended on (MM/DD/YYYY) I hereby state that I have reviewed and understand the contents of the above identified specification, including the claims, as amended by any amendment specifically referred to above. I acknowledge the duty to disclose information which is material to patentability as defined in 37 CFR 1.56, including for continuation-in-part applications, material information which became available between the filing date of the prior application and the national or PCT international filing date of the continuation-in-part application. I hereby claim foreign priority benefits under 35 U.S.C. 119(a)-(d) or 365(b) of any foreign application(s) for patent or inventor's certificate, or 365(a) of any PCT international application which designated at least one country other than the United States of America, listed below and have also identified below, by checking the box, any foreign application for patent or inventor's certificate, or any PCT international application having a filing date before that of the application on which priority is claimed.				
Prior Foreign Application Number(s)	Country	Foreign Filing Date (MM/DD/YYYY)	Priority Not Claimed	Certified Copy Attached? YES NO
			☐ ☐ ☐ ☐	☐ ☐ ☐ ☐ ☐ ☐ ☐ ☐
☐ Additional foreign application numbers are listed on a supplemental priority data sheet PTO/SB/02B attached hereto.				

91

DECLARATION - Utility or Design Patent Application		
I hereby claim the benefit under 35 U.S.C. 119(e) of any United States provisional application(s) listed below.		
Application Number(s)	Filing Date (MM/DD/YYYY)	☐ Additional provisional application numbers are listed on a supplemental priority data sheet PTO/SB/028 attached hereto.
I hereby claim the benefit under Title 35, United States Code, §120 of any United States application(s) listed below and, insofar as the subject matter of each of the claims of this application is not disclosed in the prior United States application in the manner provided by the first paragraph of Title 35, United States Code, §112, I acknowledge the duty to disclose material information as defined in Title 37, Code of Federal Regulations, §1.58(a) which occurred between the filing date of the prior application and the national or PCT international filing date of this application:		
Application Serial No.	Filing Date	Status
		Pending Pending Pending
I hereby appoint: ☒ Practitioners at Customer Number 000027777 -- AND ☐ Practitioner(s) named below: <u>Name</u> <u>Registration Number</u> Place Customer Number Bar Code Label Here		
as my/our attorney(s) or agent(s) to prosecute the application identified above, and to transact all business in the United States Patent and Trademark Office connected therewith.		
Address all telephone calls to Andrew C. Farmer at telephone number (732) 524-2925.		
Direct all correspondence to: Customer Number ☒ or Bar Code Label 000027777 OR ☐ Correspondence address below		
Name:		
Address:		
Address:		
City:	State:	ZIP
Country	Telephone:	Fax:

40
92

I hereby declare that all statements made herein of my own knowledge are true and that all statements made on information and belief are believed to be true; and further that these statements were made with the knowledge that willful false statements and the like so made are punishable by fine or imprisonment, or both, under 18 U.S.C. 1001 and that such willful false statements may jeopardize the validity of the application or any patent issued thereon.			
NAME OF SOLE OR FIRST INVENTOR:		☐ A petition has been filed for this unsigned inventor	
Given Name (first and middle (if any)) Ben		Family Name or Surname Fryer	
Inventor's Signature		Date	
Residence: City Lake Forest	State CA	Country USA	Citizenship USA
Mailing Address 25881 View Points			
City Lake Forest	State CA	ZIP 92630	Country USA
I hereby declare that all statements made herein of my own knowledge are true and that all statements made on information and belief are believed to be true; and further that these statements were made with the knowledge that willful false statements and the like so made are punishable by fine or imprisonment, or both, under 18 U.S.C. 1001 and that such willful false statements may jeopardize the validity of the application or any patent issued thereon.			
NAME OF SECOND INVENTOR:		☐ A petition has been filed for this unsigned inventor	
Given Name (first and middle (if any)) Peter C.		Family Name or Surname Zhu	
Inventor's Signature		Date	
Residence: City Irvine	State CA	Country USA	Citizenship USA
Mailing Address 8 Cresthaven			
City Irvine	State CA	ZIP 92604	Country USA

41
93

Figure 1

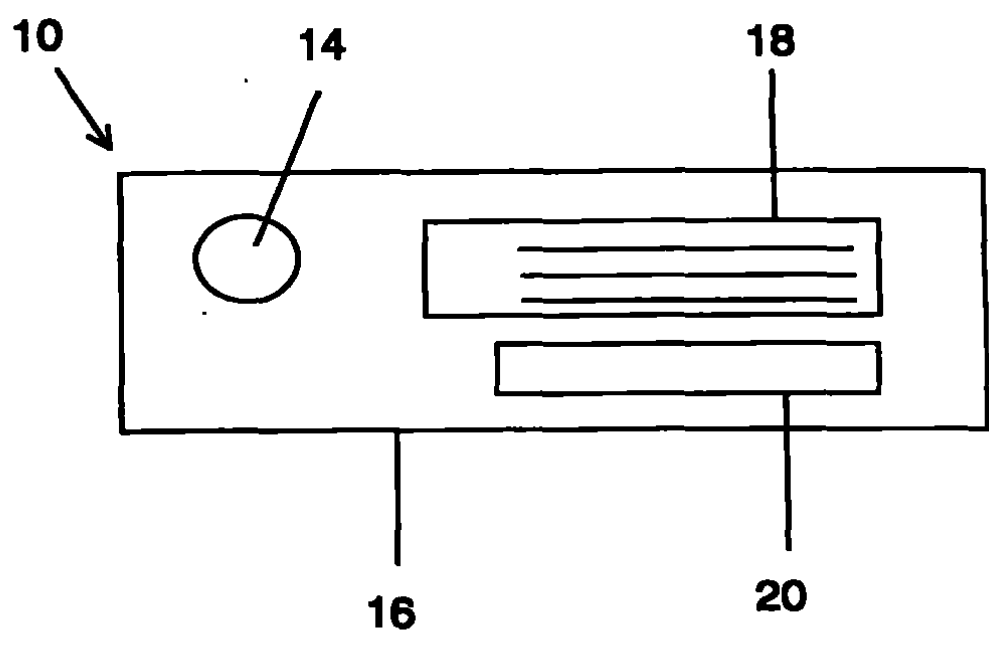
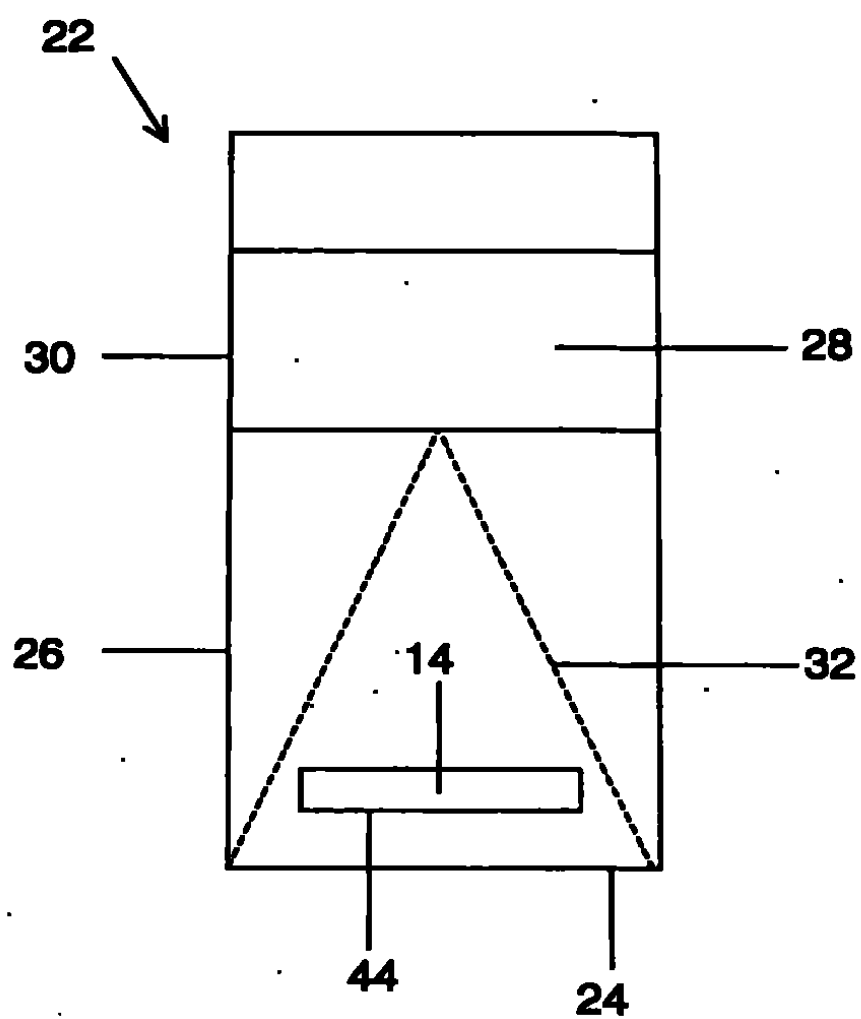


Figure 2

43
95

Figure 3

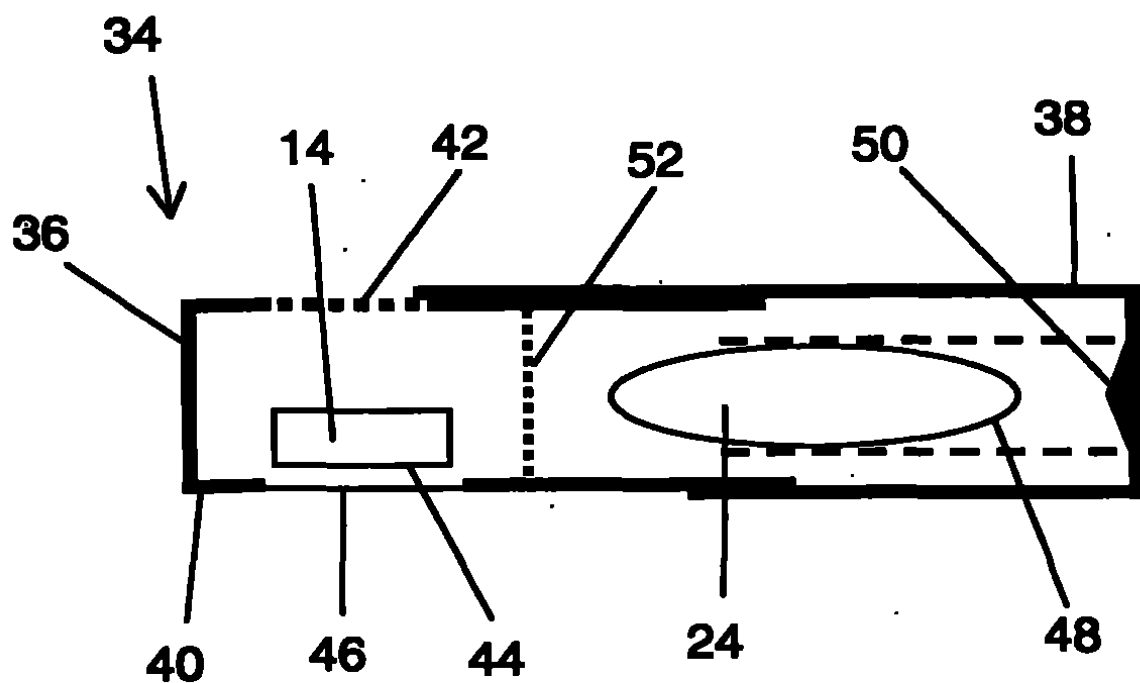


Figure 4